VE-Cadherin Rabbit PolymAb®

Catalog No.: A0734PM



Basic Information

Observed MW

90kDa,120-140kDa

Calculated MW

88kDa

Category

Primary antibody

Applications

WB,IHC-P,ELISA

Cross-Reactivity

Human, Mouse, Rat

Background

This gene encodes a classical cadherin of the cadherin superfamily. The encoded preproprotein is proteolytically processed to generate the mature glycoprotein. This calcium-dependent cell-cell adhesion molecule is comprised of five extracellular cadherin repeats, a transmembrane region and a highly conserved cytoplasmic tail. Functioning as a classical cadherin by imparting to cells the ability to adhere in a homophilic manner, this protein plays a role in endothelial adherens junction assembly and maintenance. This gene is located in a gene cluster in a region on the long arm of chromosome 16 that is involved in loss of heterozygosity events in breast and prostate cancer.

Recommended Dilutions

WB 1:2000 - 1:4000

IHC-P 1:500 - 1:2000

ELISA Recommended starting concentration is 1 μ g/mL.

Please optimize the concentration based on your specific assay requirements.

Immunogen Information

 Gene ID
 Swiss Prot

 1003/12562
 P33151/P55284

Immunogen

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

Synonyms

7B4; CD144

Contact

<u>a</u>	400-999-6126
\bowtie	cn.market@abclonal.com.cn
<u>~</u>	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS containing 50% glycerol and 0.05% BSA, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.

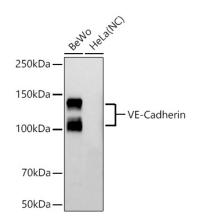

Western blot analysis of various lysates using VE-Cadherin Rabbit PolymAb® (A0734PM)at 1:4000 dilution incubated overnight at 4°C.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25 µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Negative control (NC): HeLa Exposure time: 90s.

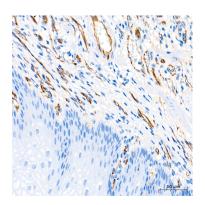


Immunohistochemistry analysis of paraffinembedded Human esophagus tissue using VE-Cadherin Rabbit PolymAb® (A0734PM) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

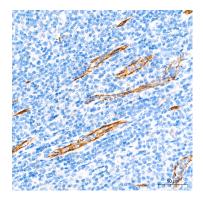


Immunohistochemistry analysis of paraffinembedded Human tonsil tissue using VE-Cadherin Rabbit PolymAb® (A0734PM) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

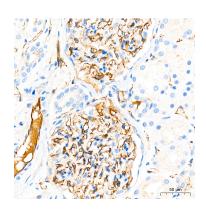


Immunohistochemistry analysis of paraffinembedded Human kidney tissue using VE-Cadherin Rabbit PolymAb® (A0734PM) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

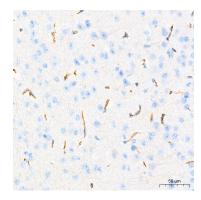


Immunohistochemistry analysis of paraffinembedded Mouse brain tissue using VE-Cadherin Rabbit PolymAb® (A0734PM) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

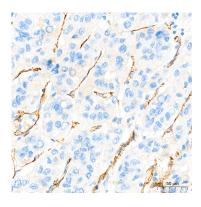


Immunohistochemistry analysis of paraffinembedded Human liver cancer tissue using VE-Cadherin Rabbit PolymAb® (A0734PM) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

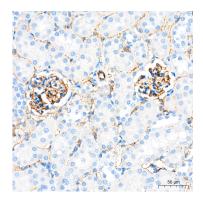


Immunohistochemistry analysis of paraffinembedded Mouse kidney tissue using VE-Cadherin Rabbit PolymAb® (A0734PM) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

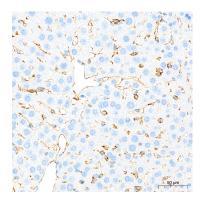


Immunohistochemistry analysis of paraffinembedded Mouse liver tissue using VE-Cadherin Rabbit PolymAb® (A0734PM) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

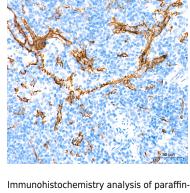


Immunohistochemistry analysis of paraffinembedded Mouse spleen tissue using VE-Cadherin Rabbit PolymAb® (A0734PM) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

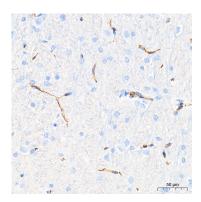


Immunohistochemistry analysis of paraffinembedded Rat brain tissue using VE-Cadherin Rabbit PolymAb® (A0734PM) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

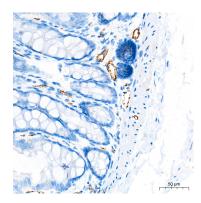

Immunohistochemistry analysis of paraffinembedded Rat colon tissue using VE-Cadherin Rabbit PolymAb® (A0734PM) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

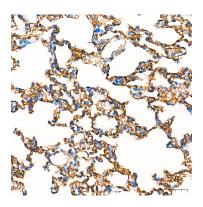

Immunohistochemistry analysis of paraffinembedded Rat lung tissue using VE-Cadherin Rabbit PolymAb® (A0734PM) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

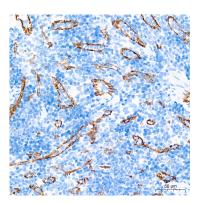

Immunohistochemistry analysis of paraffinembedded Rat spleen tissue using VE-Cadherin Rabbit PolymAb® (A0734PM) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.