LAMP3 Rabbit pAb

Catalog No.: A2895 2 Publications



Basic Information

Observed MW

70kDa

Calculated MW

44kDa

Category

Primary antibody

Applications

WB,IHC-P,ELISA

Cross-Reactivity

Human, Mouse, Rat

Background

Dendritic cells (DCs) are the most potent antigen-presenting cells. Immature DCs efficiently capture antigens and differentiate into interdigitating dendritic cells (IDCs) in lymphoid tissues that induce primary T-cell responses (summary by de Saint-Vis et al., 1998 [PubMed 9768752]).

Recommended Dilutions

WB 1:1000 - 1:5000

IHC-P 1:50 - 1:100

ELISA Recommended starting concentration is 1 μ g/mL. Please optimize the

concentration based on your specific assay requirements.

Immunogen Information

Gene ID27074

Swiss Prot
Q9UQV4

Immunogen

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

Synonyms

LAMP; CD208; DCLAMP; LAMP-3; TSC403; DC LAMP; DC-LAMP; LAMP3

Contact

<u>a</u>		400-999-6126
\bowtie		cn.market@abclonal.com.cn
•	T	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS containing 50% glycerol, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.

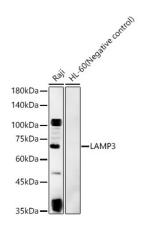


Western blot analysis of various lysates, using LAMP3 Rabbit pAb (A2895) at 1:3000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (RM00020).

Exposure time: 180s.

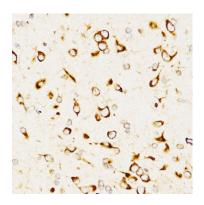

Immunohistochemistry analysis of paraffinembedded Rat brain using LAMP3 Rabbit pAb (A2895) at dilution of 1:200 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

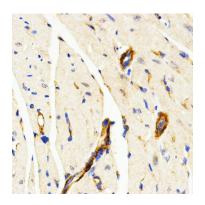

Immunohistochemistry analysis of paraffinembedded Mouse heart using LAMP3 Rabbit pAb (A2895) at dilution of 1:200 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.