UBIAD1 Rabbit pAb

Catalog No.: A7640



Basic Information

Observed MW

37kDa

Calculated MW

37kDa

Category

Primary antibody

Applications

WB,IHC-P,ELISA

Cross-Reactivity

Human, Mouse, Rat

Background

This gene encodes a protein thought to be involved in cholesterol and phospholipid metabolism. Mutations in this gene are associated with Schnyder crystalline corneal dystrophy.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:100 - 1:200

ELISA Recommended starting concentration is 1 μg/mL.

Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID29914

Swiss Prot
Q9Y5Z9

Immunogen

Synthetic peptide. This information is considered to be commercially sensitive.

Synonyms

SCCD; TERE1; UBIAD1

Contact

a		400-999-6126
\bowtie		cn.market@abclonal.com.cn
\overline{a}	ı	www.ahclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification


Storage

Store at -20°C. Avoid freeze / thaw cycles.


Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

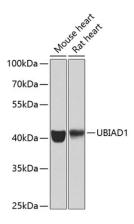


Western blot analysis of various lysates using UBIAD1 Rabbit pAb (A7640) at 1:1000 dilution._Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution._Lysates/proteins: 25µg per lane._Blocking buffer: 3% nonfat dry milk in TBST._Detection: ECL Enhanced Kit (RM00021)._Exposure time: 60s.

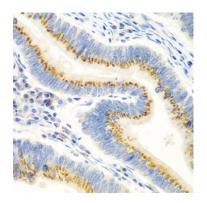

Immunohistochemistry analysis of paraffinembedded Human colon carcinoma using UBIAD1 Rabbit pAb (A7640) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

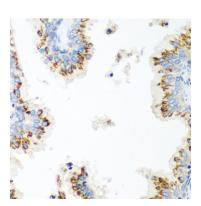

Immunohistochemistry analysis of paraffinembedded Human prostate using UBIAD1 Rabbit pAb (A7640) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

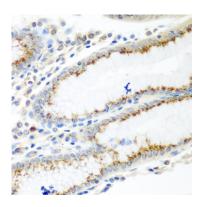

Immunohistochemistry analysis of paraffinembedded Human stomach using UBIAD1 Rabbit pAb (A7640) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.